




FAMED SAĞLIK TIBBİ CİHAZLAR
İTH. İHRC. SAN. ve TİC. LTD. ŞTİ.

HD-1315P Pro

New generation
all-in-one
portable & multi-
functional
endoscope system

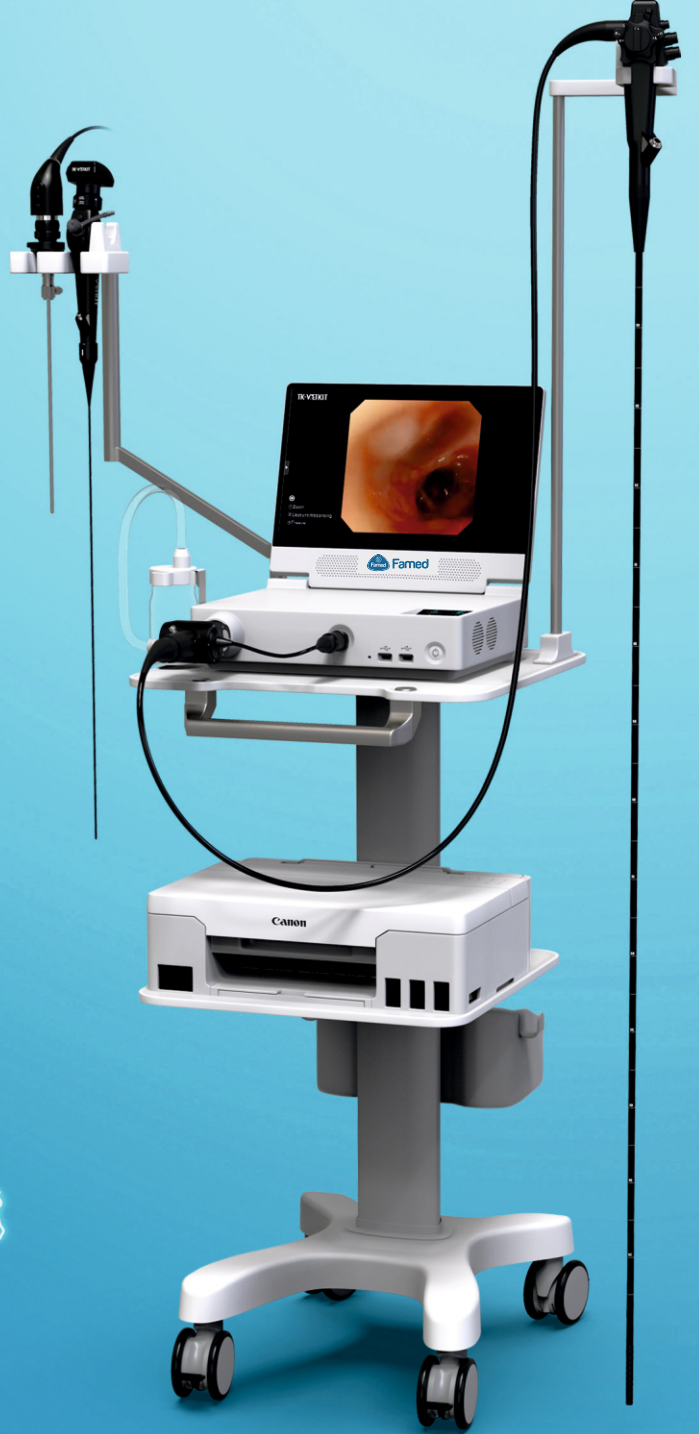


www.famed.com.tr



info@famed.com.tr

Full HD (1080P) image processor for for high-quality image visualization

Experience the cutting-edge technology of FAMED Full HD high-definition CMOS sensor, delivering unparalleled image quality and an optimized field of vision. With this advanced technology, veterinarians can now see with remarkable clarity and diagnose with exceptional precision during endoscopic surgeries.,

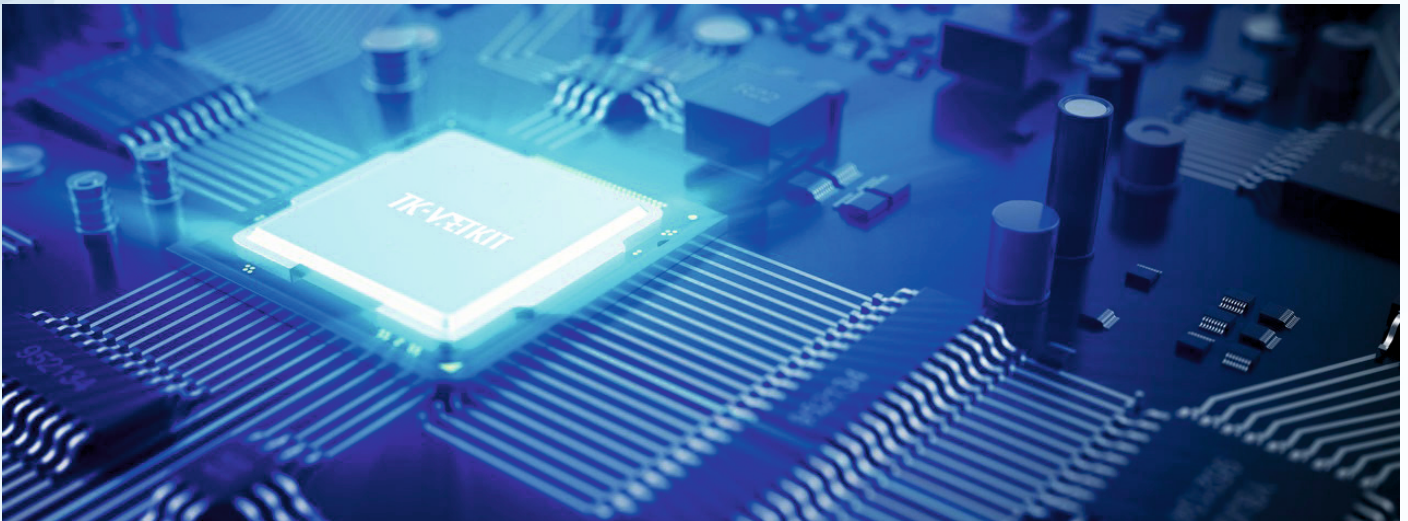

Portable and foldable design with removable batteries for both indoor and outdoor scenes



All in one foldable design



Inside battery last for 1.5 h



Lighter weight,
easy to carry

Multi-scenario applications



Installed in hospital



Suitable for outdoor



Easy for transportation

Compatible with rigid endoscopes and flexible endoscopes

The HD-1315P Pro is the perfect companion for FAMED endoscopes from the EU Series, GU Series, and RS Series. This cutting-edge device is designed to be compatible with a wide range of veterinary surgical and diagnostic applications, ensuring versatility and precision in every procedure



EU Series



- Bronchoscopy
- Rhinoscopy
- Vaginoscopy
- Otoscopy
- Cystoscopy
- Small animals gastroscopy

Model	Resolution	Tube	Working Channel	Working Length	Up/Down	L/R Rotation
FMD-120	720*720	2.8mm	1.2mm	610mm	U210°/D130°	L120°/R120°
FMD-130	720*720	3.8mm	1.5mm	1000mm	U180°/D130°	L120°/R120°
FMD-140H	HD 720P	4.8mm	2.0mm	1000mm	U180°/D130°	L120°/R120°
FMD-250H	HD 720P	5.8mm	2.8mm	1500mm	U180°/D130°	L120°/R120°

GU Series



- Gastrointestinal endoscopy
- Duodenoscopy
- Esophageal endoscopy
- Colonoscopy
- Gastroscopy

Model	Outer Diameter	Working Channel	Working Length	Up/Down	L/R Rotation	View Angle	Function Keys	Water Jet	Resolution
FMD-160H	6.8mm	2.2mm	1500mm	U210°/D130°	L100°/R100°	120°	1	/	HD
FMD-170Q	7.8mm	2.8mm	1500mm	U210°/D130°	L100°/R100°	120°	3	/	HD
FMD-170FH	7.8mm	2.2mm	1500mm	U210°/D130°	L100°/R100°	140°	3	/	FHD
FMD-190FQ	9.2mm	2.8mm	1500mm	U210°/D130°	L100°/R100°	120°	1	/	FHD
FMD-190FHPL	9.2mm	2.8mm	3000mm	U180°/D130°	L100°/R100°	140°	3	Yes	FHD
FMD-1B0FHPL	11.5mm	3.2mm	3500mm	U180°/D130°	L100°/R100°	140°	3	Yes	FHD



- Laparoscopy
- Arthroscopy
- Otoscopy and Rhino-laryngoscopy
- Cystoscopy

Item	Model	Description
Camera Head	CM-U12	Full HD 1080P

15.6" FullHD monitor with Medical Record Management System for better performance

- Utility functions: Easily capture photos and take videos and view important details by freezing images.
- Medical report: Easily access features, patient information, user settings, and system configuration through the management system.
- Customizable buttons – Customize endoscope buttons for functions like taking photos, taking videos, white balance, freezing and adjusting brightness.



Multi levels LED light source and 3 levels air/water pump



Meet various needs for animals in different sizes to do air inflation and adjust suitable brightness in the cavity

Pump Level 1



Pump Level 2



Pump Level 3





 | www.famed.com.tr



FAMED SAĞLIK TIBBİ CİHAZLAR
İTH. İHRC. SAN. ve TİC. LTD. ŞTİ.



Address:: Cevizli Mh. Hacılar Cd. No:44/A Kartal - İstanbul / TURKEY
Phone : +90 216 597 0 567 / +90 216 441 29 99 **Fax:** +90 216 397 15 15
W: www.famed.com.tr **E:** info@famed.com.tr